


Give every newborn the best start in life



## Ecosy630 Infant Warmer

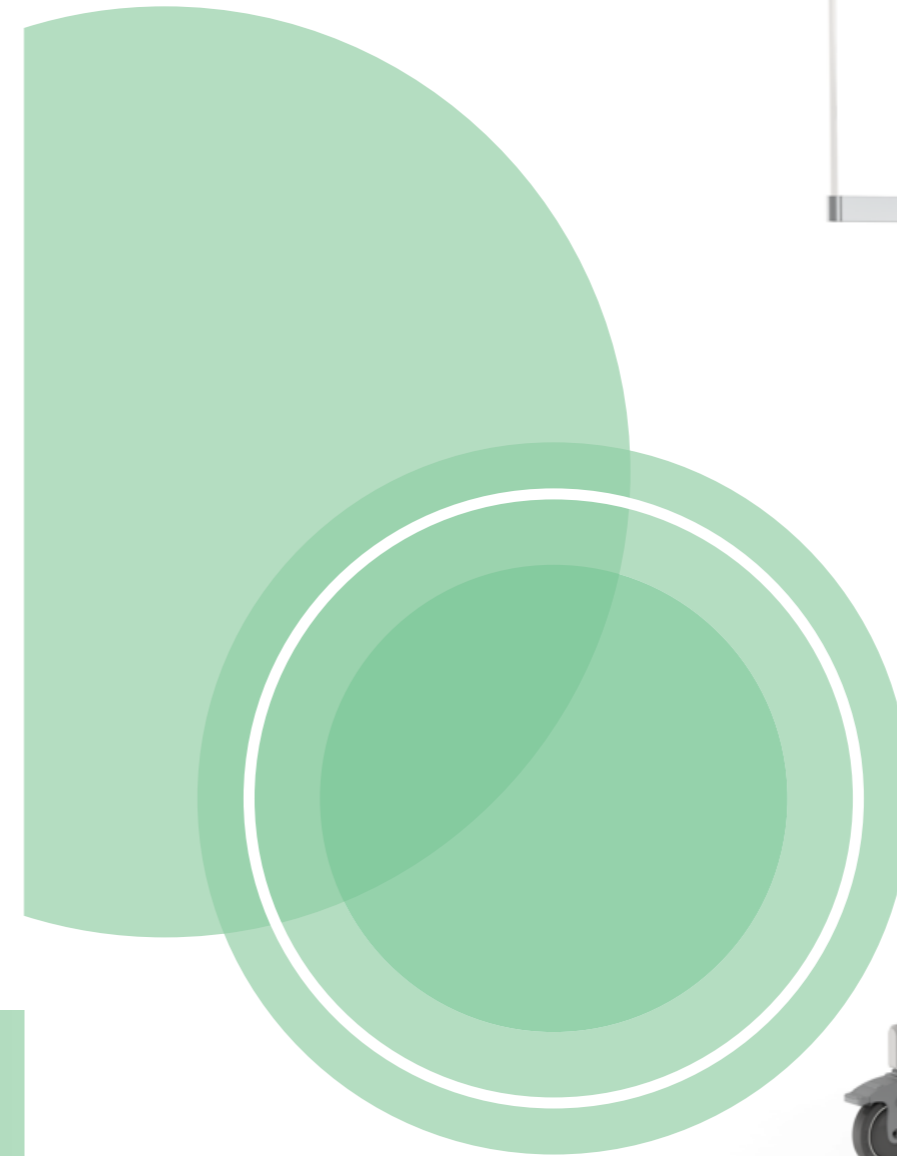

### BNG MEDICAL INSTRUMENTS CO.,LTD.

📍 3rd Floor, Block 3, No.14, Zhongxing Road, Xiuxin Community,  
Kengzi Street, Pingshan District, Shenzhen, China  
✉ Email: [info@bngmed.com](mailto:info@bngmed.com)      ☎ Tel: +86-0755-27392867  
🌐 Web: [www.bngmed.com](http://www.bngmed.com)


CE ISO  
2460 13485

[www.bngmed.com](http://www.bngmed.com)

# Ecosy 630 INFANT WARMER

## Basic configurations

- Far infrared quartz heating tube
- Baby, manual and pre-heating modes
- APGAR and CPR timer
- Transparent guard door on four sides
- Baby bed tilting angle from 0~13°

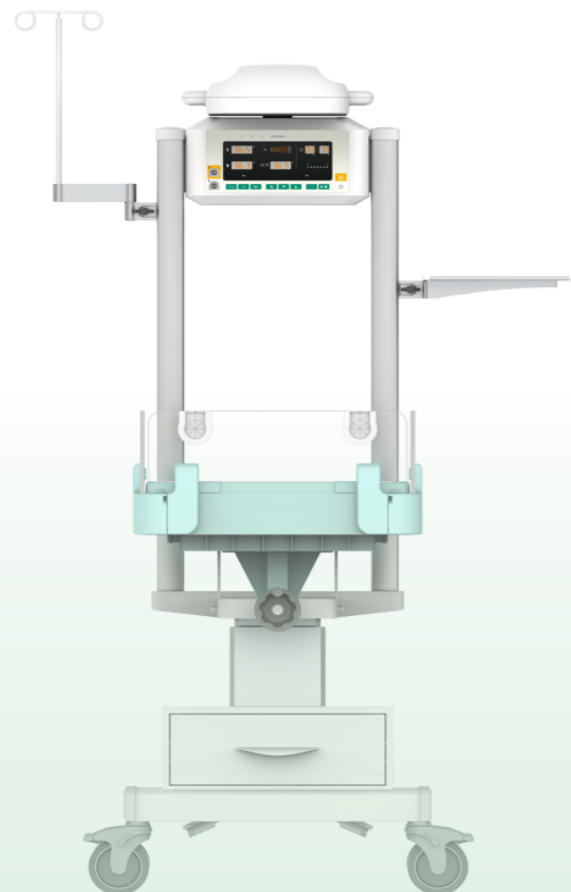


## Optional configurations

- Storage drawer
- X-ray cassette
- Lifting stand


## ■ Heating head module

Heating head adjustable angles:  $\pm 120^\circ$  horizontal rotation,  $\pm 30^\circ$  vertical rotation.

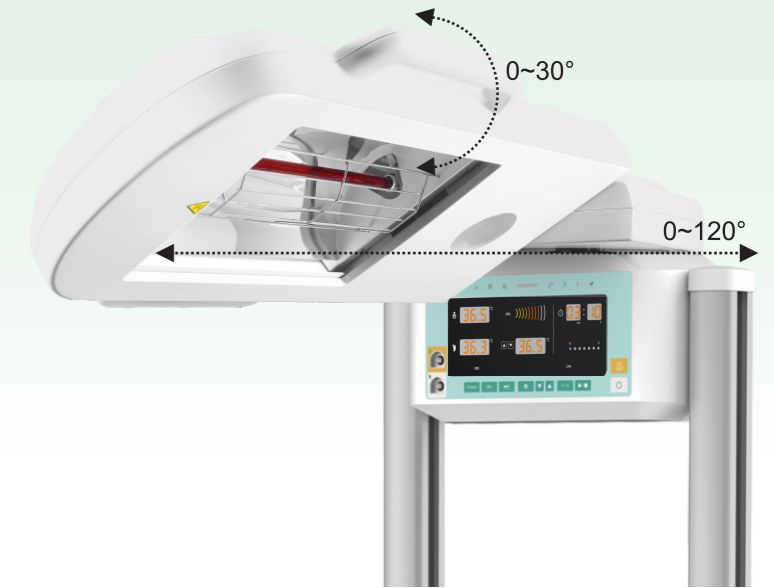

## ■ LED display

Easily to operate the device and check information.


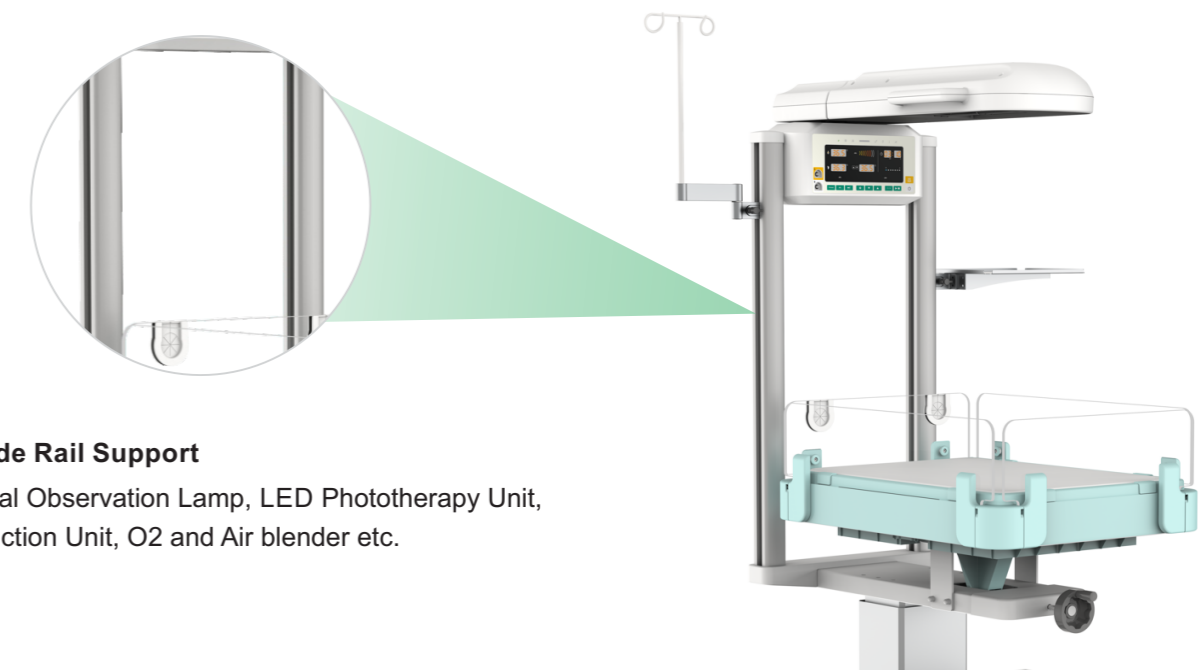


## ■ Multi-functional Side Rail Support

Design for the external Observation Lamp, LED Phototherapy Unit, Resuscitator Unit, Suction Unit, O2 and Air blender etc.

■ **Optional storage drawer**

Storage drawer equipped with damping technology for silent operation and can be operated on both left and right sides which is ideal for storage.

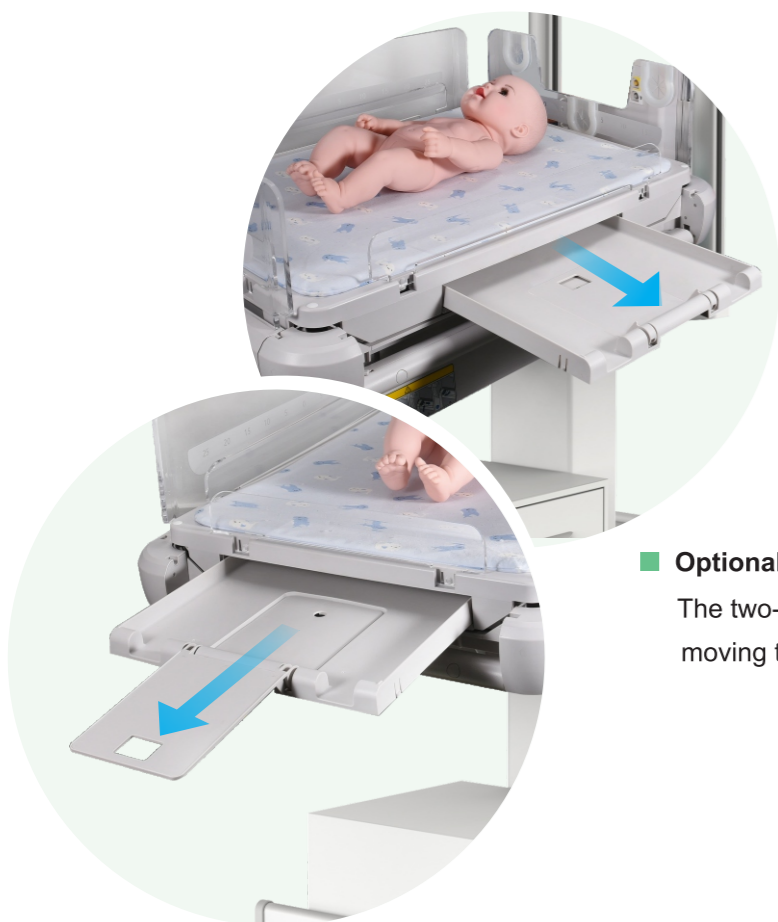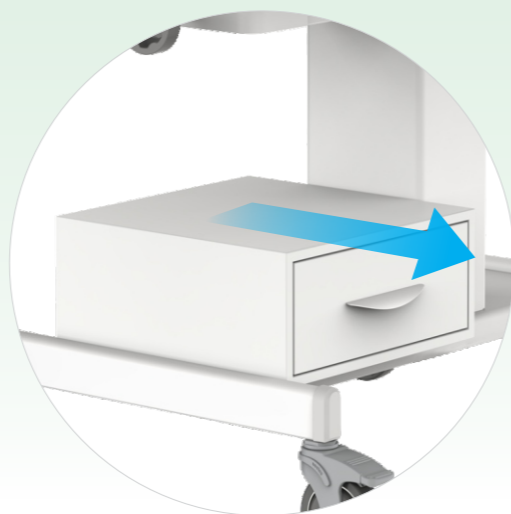

■ **Optional X-Ray cassette**

The two-position pull-out X-ray cassette under the baby tray, avoid moving the newborns when taking the X-Ray examination.

■ **Lifting stand**

Adjust the height of the infant radiant warmer based on operator's height. The foot pedal can be operated from both left and right sides.



## Ecosy 630 Infant Warmer

Heater capacity	500W
Temperature control	Servo(Baby) Mode/Manual/Preheating mode
Skin temperature	Setting range:34.0~38.0°C
	Display range:20.0~42.0°C
Heater output	Setting range:0~100%(in 5% increments)
	Indication range:0~100%(in 1% increments)
Timer	APGAR mode/CPR mode switchable
	APGAR Beep:1,5,10 min
Display	LED digital display
Bed tray	Inclination angle(back and forth):±13°
	Mattress surface:500x700x25mm
Dimension	840(L)*1100(W)*1790(H)mm

## Optional Functions

X-ray cassette	
Dimension	53x35cm
Height adjustment	Up to 20cm
Storage/Shipping Information	
Temperature	-25~55°C
Humidity	0~93% RH

Notice: Specifications subject to changes without prior notice. All rights reserved by BNG Medical.

